

Activate the Sensopro integration

1. Make sure integrations are included within your Sensopro subscription
2. Activate the Sensopro integration in the Magicline Appstore
3. Enter the desired Sensopro Luna subscriptions or admissions in Magicline
4. Make sure you have accepted the data exchange regulations
5. Make sure the Sensopro Luna has an internet connection
6. Sensopro will contact you to schedule the activation
7. Sensopro activates the integration on your Sensopro Luna

Set up Sensopro subscriptions and entry passes

The screenshot shows the Magicline app interface. The top navigation bar includes the Magicline logo, a search icon, and the text 'OpenApi' and 'Sensopro'. The main content area is divided into three tabs: 'Einstellungen', 'Integrationen', and 'Partnerintegrations'. The 'Partnerintegrations' tab is selected. Below the tabs is a search bar and a list of integration categories. The 'Partnerintegrations' category is highlighted in blue. Below this, there is a table with two columns: 'NAME' and 'STATUS'. The table contains one entry: 'Sensopro' with a status of 'Aktiv'. A red arrow points to the 'Aktiv' status.

NAME	STATUS
Sensopro	Aktiv

1. Once the Sensopro integration has been activated in the app store it shows up under:

Settings -> Integrations -> Partner integrations

Set up Sensopro subscriptions and entry passes

2. Under «Add inclusive contingent», add the Sensopro benefit to your desired subscriptions.

The screenshot shows the Magicline management interface. The left sidebar contains a navigation menu with items like Dashboard, Check-In, Mitglieder, Interessenten, Finanzen, Termine, Kurse, Kommunikation, Analytics, Verkauf, Aufgaben, Ressourcen, Trainingswelt, Einstellungen, Hilfe und Support, Einkaufen, and Appstore. The main content area is titled 'Partnerintegrationen' and shows a configuration for 'Sensopro Aktiv'. The 'Produkte' tab is selected, and the 'sensopro' entry is visible. A red arrow points to the three-dot menu icon next to the 'sensopro' entry. Below the 'sensopro' entry, there are two buttons: 'INKLUSIVKONTINGENT HINZUFÜGEN' and 'KARTENANGEBOTE HINZUFÜGEN'. A second red arrow points to the 'INKLUSIVKONTINGENT HINZUFÜGEN' button.

The screenshot shows the 'Inklusivkontingent hinzufügen' dialog box. The dialog has a title bar with a close button (X). The main content area contains several dropdown menus: 'Volumen' (Flat), 'Nutzungszeitraum' (Einmalige Nutzung), and 'Gilt für' (Tarife). Below these, there is a dropdown menu with the text 'Bitte wählen...' and a list of options: 'Basic', 'Basic', 'Premium', and 'Premium Abo'. A red arrow points to the 'INKLUSIVKONTINGENT HINZUFÜGEN' button from the previous screenshot.

Set up Sensopro subscriptions and entry passes

The screenshot shows the Magicline Sensopro configuration interface. The left sidebar contains a navigation menu with items like Dashboard, Check-In, Mitglieder, Interessenten, Finanzen, Termine, Kurse, Kommunikation, Analytics, Verkauf, Aufgaben, Ressourcen, Trainingswelt, Einstellungen, and Hilfe und Support. The main content area is titled 'Sensopro Aktiv' and shows configuration options for 'Produkte' and 'Konfiguration'. Under 'Produkte', there is a section for 'Inklusive Kontingent' with a 'Premium Abo' and a 'Flat' option for 'Einmalige Nutzung'. Below this, there is a section for 'Kartenangebote' with two blue links: 'INKLUSIVKONTINGENT HINZUFÜGEN' and 'KARTENANGEBOTE HINZUFÜGEN'. A red arrow points to the 'KARTENANGEBOTE HINZUFÜGEN' link.

3. Under «add entry pass», add your desired entry passes.

The screenshot shows the 'Kartenangebote hinzufügen' dialog box. It contains the following fields and options:

- Nutzungsumfang:** 1
- Preis:** 0,00 €
- Nutzungszeitraum:** Einmalige Nutzung
- Artikelnummer:** (empty field)
- Gültigkeit:** Unbegrenzt gültig

At the bottom right, there are two buttons: 'ABBRECHEN' and 'HINZUFÜGEN'.


Set up Sensopro subscriptions and entry passes

The screenshot displays the MAGICLINE interface for configuring Sensopro. The left sidebar contains navigation options: Dashboard, Check-In, Mitglieder, Interessenten, Finanzen, Termine, Kurse, Kommunikation, Analytics, Verkauf, Aufgaben, Ressourcen, Trainingswelt, Einstellungen, Hilfe und Support, and Einkaufen. The main content area is titled 'Sensopro' and shows configuration options for 'Produkte' and 'Konfiguration'. Under 'Produkte', there are two sections: 'Inklusive Kontingent' and 'Kartenangebote'. The 'Inklusive Kontingent' section lists 'Premium Abo' with 'Flat' and 'Einmalige Nutzung' options. The 'Kartenangebote' section lists two types of passes: '1 x sensopro Art-Nr. V-3 Einmalige Nutzung' and '10 x sensopro Art-Nr. V-4 Einmalige Nutzung', both priced at 0,00 € and offering 'Unbegrenzte Gültigkeit'. Below these sections are buttons for 'INKLUSIVKONTINGENT HINZUFÜGEN' and 'KARTENANGEBOTE HINZUFÜGEN'. Red arrows point from the text labels on the right to the corresponding product rows in the interface.

Example subscription

Example entry pass

Log-in on the device

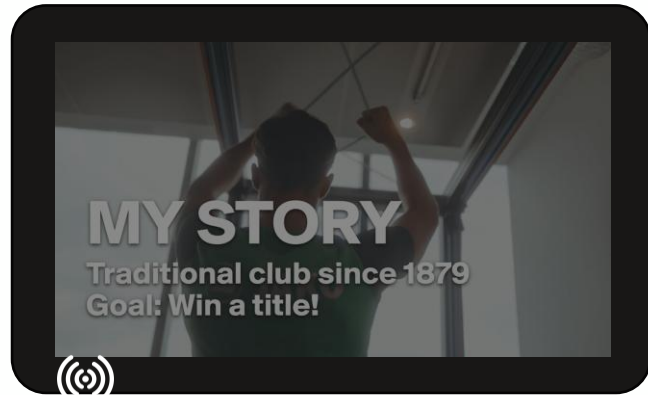

Tap the screen with
your RFID chip card



Screen unlocks:
Valid subscription / entry pass



Screen doesn't unlock:
Invalid subscription / entry pass

Further informations

Accept data exchange

Wir haben unseren Auftragsverarbeitungsvertrag gemäß Art. 28 DSGVO für Magicline aktualisiert. Bitte bestätige diese, damit du Magicline weiterhin nutzen kannst. [Lesen & akzeptieren](#)

MAGICLINE
Suche
OpenApp Admin User

Einstellungen Integrationen Aktivierte Integrationen Sensopro

Suche Sensopro aktiv
Aktivierungsdatum: 09.02.2024
Entstehungsdatum: 09.02.2024

Produkte Konfiguration **Systemschutz**

Berechtigungen

Produkt	Konfiguration	Status	Systembestandteil	Status
CUSTOMER_BENEFIT_READ		Aktiviert	System bestehender Mitglieder teilen	Aktiviert
CUSTOMER_BENEFIT_WRITE		Aktiviert	System bestehender Interessanten teilen	Aktiviert
CUSTOMER_READ		Aktiviert	System ehemaliger Mitglieder freigeben	Deaktiviert
			System neuer Mitglieder teilen	Aktiviert
			System neuer Interessanten teilen	Aktiviert

Systembestandteil mit Partnerintegrationen

Webhook Events

Use with «Sensopro Hub» to assign trainingplans and track trainings and coordination tests.

SENSOPRO

HEY PHYSIO ATHLETIK

Suche Trainingsplan

Trainingsplan
Max Muster

Max Muster
0022344566666

COORDINATION SCORE ANALYSE TRAININGPLAN TRAININGSERIE

Coordination Score
LEISTUNGS SCORE
71.5 / 100

Schwindigkeit
Rhythmus
Gleichgewicht
Stabilität

Schwindigkeit: 65.3
Gleichgewicht: 70.1
Stabilität: 88.2
Rhythmus: 88.2

71.5 / 100
75
TRAINING HINZUFÜGEN